

# X-ray table Amadeo V pico for the small animal practice

The inexpensive Amadeo V pico X-ray table was specially developed for small animal practices. It has an integrated stand including a 5 kW generator and can be used in a variety of ways. With the Amadeo V pico, you can quickly and easily take exposures with the detector placed above or under the table. The moving cassette drawer makes it easier to position the four-legged patient.



- X-ray unit can be moved longitudinally centred over the centre of the table, cassette / grid drawer follows the generator in a centred position
- Standardised with 5 kW generator, other generators on request
- Table frame is made of anodised aluminium
- Table top is made of radiolucent white plastic
- Easy to clean and resistant to disinfectants
- Optional: size-adjustable HD fixing cot made of plastic
- Optional: also available without generator
- Table top dimensions: 1,591 x 742 mm
- Working height: 806 mm; total height: 2,215 mm

# X-ray table Amadeo V pico






Generator with large touchscreen

Simple & uncomplicated operation



[other generators on request]

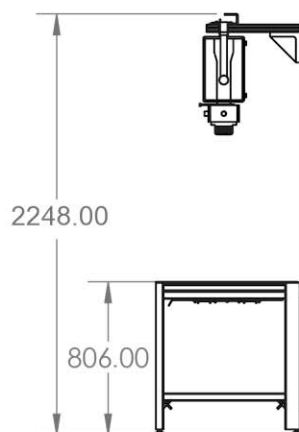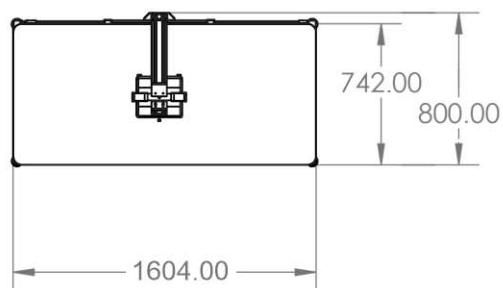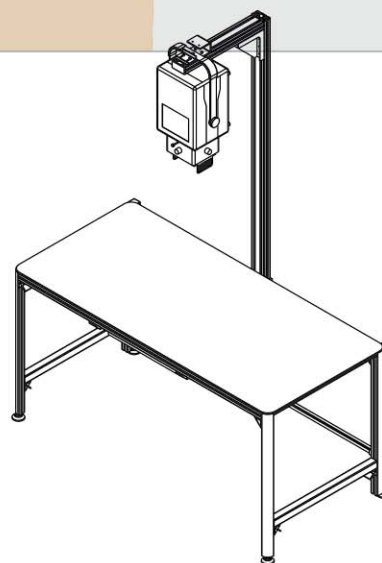

Please visit [www.or-technology.com](http://www.or-technology.com) for more details

## Specifications subject to revision without notice

The editor strives to impart correct and up to date information. The provided specifications are based on current knowledge and are subject to revision without notice. This brochure is subject to correction. The editor assumes no responsibility for the information being up to date, correct and complete. All furnished logos, pictures and graphics are property of the particular company and subject to copyright of the licensor. Use, dissemination, distribution or copying of the pictures, logos or text compiled or processed by the editor is subject to our written consent. All rights reserved.

**Headquarters:** OR Technology (Oehm und Rehbein GmbH), 18057 Rostock, Germany, Neptunallee 7c  
Info-Hotline: +49 381 36 600 600, [info@or-technology.com](mailto:info@or-technology.com), [www.or-technology.com](http://www.or-technology.com)  
**OR Technology UK:** Celtic SMR Ltd., Frederick House, Hayston View, Johnston, Haverfordwest, Pembrokeshire SA62 3AQ, United Kingdom, [www.celticsmr.co.uk](http://www.celticsmr.co.uk), [sales@celticsmr.co.uk](mailto:sales@celticsmr.co.uk)